


chromotek
new tools for better research

NANO-TRAPS

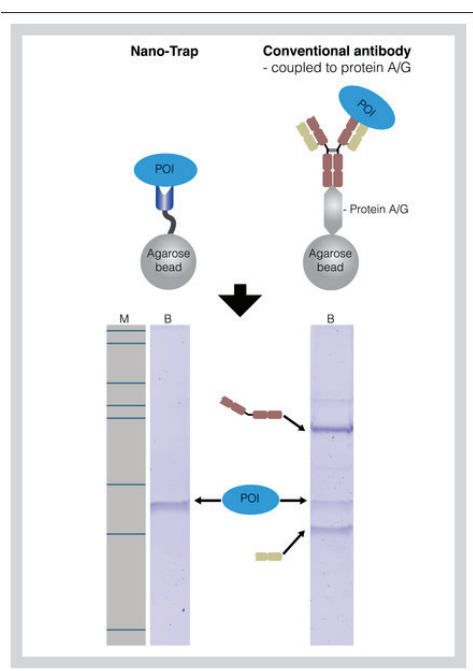

Nano-Traps are ideal for **fast, reliable and efficient one-step immunoprecipitations.**

They pull down your proteins of interest (POI) and their interacting factors from cell extracts or organelles.

Nano-Traps consist of a VHH, also termed nanobody, coupled to an immobilizing matrix (see Technology).

You can use Nano-Traps for a multitude of biochemical analyses such as:

- Immunoprecipitation (IP) / Co-IP
- Mass spectrometry
- Enzyme activity measurements
- CHIP / RIP analysis

Benefits

- No heavy and light chains in your downstream gel and mass spectrometry analysis
- High binding affinity: Dissociation constants down to 1pM
- Short incubation times; reliable and consistent results
- More than 1.000 publications
- Fulfills highest requirements on antibody validation: Structure and function characterized



Case rappresentate

abcam

Abnova

ACEA
Biosciences, Inc.

Biomedal

BIO **ROTH**

Caprico
Biotechnologies

chromotek
new tools for better research

eBioscience
ThermoFisher
SCIENTIFIC

GeneTex

Goodfellow

IBL **KPL**

LEXOGEN

MERCK MILLIPORE **CHEMICON**
MERCK **upstate**
cell signaling solutions

Percipio **oxis**
Biosciences

RayBiotech
Empowering your proteomics

SERVA
Electrophoresis

SB Sino Biological
210 03 001 5011 001 5011 001 5011

SPECTRUM
LABS.COM

Technology

Beside conventional antibodies, Camelidae (alpacas, llamas, camels and dromedaries) possess a second type of antibodies called heavy chain antibodies (hcAbs). HcAbs are devoid of light chains and bind their antigen via a single variable domain (VHH), also known as nanobody. These VHH domains have excellent binding properties and can be produced at constant high quality without batch-to-batch variations.



Nano-Traps compared to conventional antibodies at a glance:

Nano-Traps	Conventional antibodies
Ready to use	Additional incubation step for binding of antibodies to protein A/G beads
No heavy & light antibody chains in your downstream application	Contaminating heavy & light antibody chains may appear in your downstream application
Very high binding affinities	Low to high binding affinities
Harsh wash conditions may be applied for selective protein-protein interactions	Limited buffer/reagent compatibility of IgG
Short incubation times of 5-30 minutes	Considerable incubation times of twice minimum 1 hour; this may interfere with low affinity interaction partners if you perform Co-IP
High reproducibility due to recombinant expression with very low batch to batch variation	Reproducibility can be high dependent on monoclonal – but polyclonal can be different
Validated and well characterized	Validation status depends on individual antibody
In general not suitable as detection antibody in Western Blot; however there are exceptions from this lecture books' rule	Can generally be used as detection antibody in Western Blot

Scopri di più su www.chromotek.com

Hai bisogno di altre informazioni? Contattaci!

PRODOTTI GIANNI
Via Quintiliano 30
20138 Milano Italy

Num Verde: 800 832027
Fax: 02 5097276
ricerca.assistenza@prodottigianni.com

WWW.RICERCA.IT